


Mini Laminar Flow Cabinet Mini Fume Hood

◎ Table type clean equipment.

Characteristics:

- HEPA filter with efficiency: 99.995% at 0.3um.
- Microprocessor control system, LED Display.
- Manual front window is 5mm toughened glass, anti-oxidation.
- Air speed adjustable.
- Desktop type, with a handle, saving space and easy to move.
- Pre-filter is polyester fiber, washable and anti-ultraviolet radiation.
- Automatic cutting off UV lamp function by opening glass.
- CJ-600PII, CJ-600NII have color screen LCD controller (with sterilization timing, built-in socket controllable switch, front window status display, can set three position interlock, cumulative work timing, remote control)

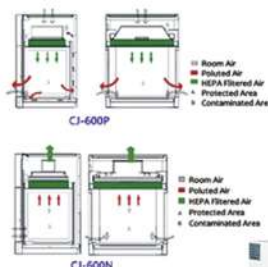

Patent No. : ZL 2018 3 0095053.5

Main Technical Parameters:

Model		CJ-600P	CJ-600N
		CJ-600PII	CJ-600NII
Applicable Scope		Protect experiment samples	Protect operator and environment
Airflow Direction		Vertical	
Function	Cleanliness	Class 100	
	High Efficiency Filter	0.3um partical ≥99.999%	
	Colony count	≤0.5/vessel (petri dish is ϕ 90mm)	
	Wind Speed	initial: 0.6;final:0.2 (m/s)	>0.5 (m/s)
	Wind Speed Range	0.2-0.6 (m/s)	
	Noise Level	≤62db	
	Vibration Semi-peak Value	≤3um	
	Work Surface Bearing	10kg	
Composition	Inner Chamber	SUS304 Stainless steel plate	
	Outer Shell	Cold rolling steel electrostatic spraying exterior	
	Operation Door	Tempered Glass	
	Operation Door Biggest Opening	310mm	
	Pre-Filter	Polyurethane	
	LED Lamp	8W*1	
	UV Lamp	8W*1	
	Top Fan	90W one piece	
Controller	Display	LED (Model II: LCD color screen)	
	Wind Speed Control	Intelligent multi-level adjustable	
	Additional Function	One key to open and close lighting sterilization	
Specification	Inner Size (W*L*H) (mm)	480*340*370	
	Exterior Size (W*L*H) (mm)	550*460*700	
	Packing Size (W*L*H) mm	680*590*850	
	Max Power	110W	
	Current Rating (50/60Hz)	AC220V/0.5A	
	NW/GW kg	36/41	
Optional Accessories		CJ-600PII, CJ-600NII: wind speed sensor	