

Shaking Water Bath

쉐이킹항온수조



Features

- The grounded circuit prevent a short circuit and over current of electricity.
- Time set function and constant temperature.
- Unique design multi rack separable. (utility model)
- Minimize noise to solve the problem of old products.
- Equipped internal bottom for the safety of the user and protection of the unit.

Specifications

Made in Korea


Spec. / Model		SWB-25	SWB-35	SWB-55
Type		Digital		
Controller		Digital PID Controller		
Range		Ambient +5 ~ 98 °C		
Temp	Accuracy	±0.2 °C		
Size (with heater cover)	Internal (W x D x H, mm)	350x300x240	500x300x240	660x350x240
	External (W x D x H, mm)	580x360x400	730x360x400	890x410x400
Safety		Over Temp limiter		
Material	Internal	Stainless steel STS #304		
	External	Steel, Powder coating		
Electric Requirements		220VAC, 50 or 60Hz		
Rack		Spring wire rack		
Heater		1500W	2300W	2800W
Cover		Stainless steel STS #304		

Shaking Water Bath


www.hysc.co.kr

쉐이킹항온수조

Distinction

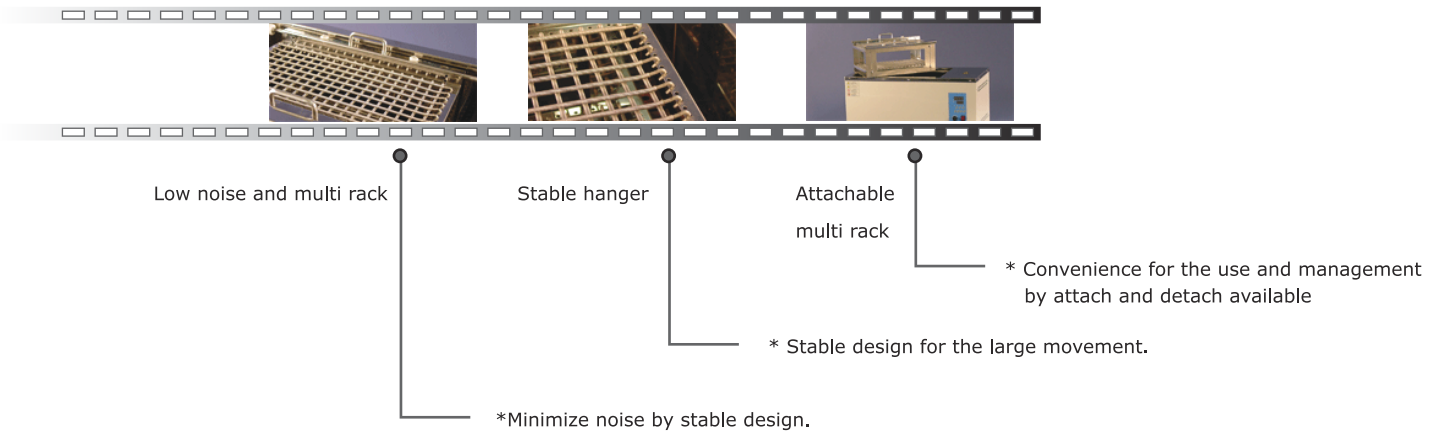

Structure

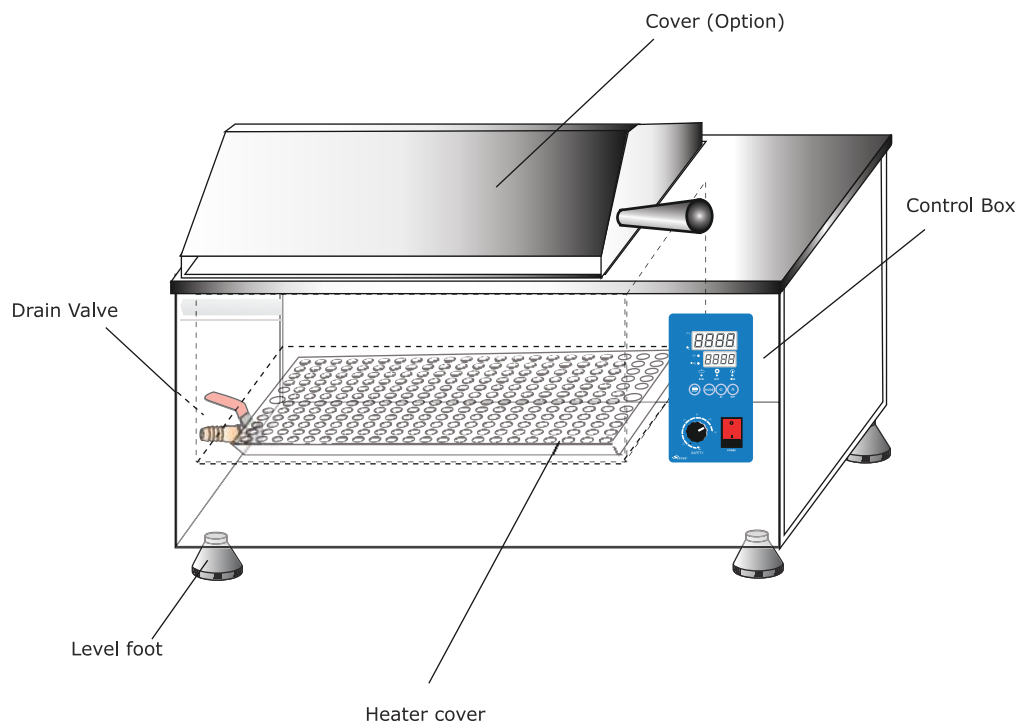

SHAKING WATER BATH

Standard equipment for a physics and chemistry laboratory

Stable temperature, horizontal reciprocating motion

Convenient removable rack



Cover


STS #304

Spring Rack / Flask holder

PID Controller

Powder coating

Over Temp Limit

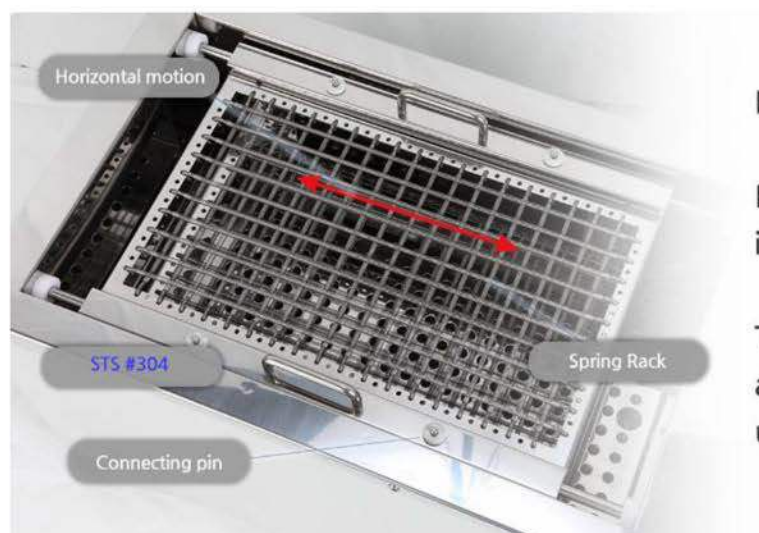

Horizontal motion

STS #304

Spring Rack

Connecting pin

Precision temperature control by PID control

Horizontal movement control
in the range of 30 - 270 rpm

The spring rack is removable and replaceable,
and can be replaced for management and
usage purposes.

