

Rotary Evaporator

HS-2005S (CE)

Distinguished Features


- 6 different types of flasks are attachable (From 50 mL to 3,000 mL)
- Digital speed control
- Digital R.P.M Display/Speed control range : 20~280 R.P.M
- Slant type condenser
- Digital interior temperature measure and display
(It catches the temperature changes caused by the interior pressure drops)
- Motorized UP-DOWN control switch
- Easy to utilize standard distillation solution
- Auto timer function
- Autonomous condenser lifting function after a complete shut down of the machine

Included Water Bath HS-3001

(PTFE Coating)

If the amount of water goes down under the certain threshold in the heating bath, the system automatically turns off its power for the safety purpose

- Included vacuum controller
- *Vacuum range : 0~760 mmHg

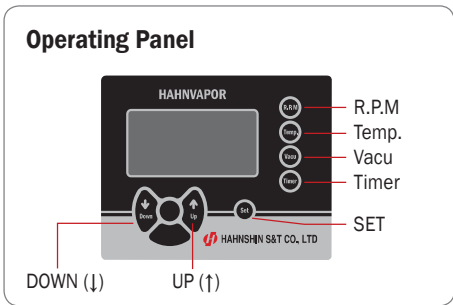

Registration No.
20-0461039/30-0668181

HS-2005S (CE)

Heater	1,000 W
Bath Material	SUS 304
Weight	2.8 Kg
Bath Dimension (Ø × H)	245 × 120 mm (Capacity : 4.3 Liter)
Required Power Source	AC 220 V, 3.5 A, 50/60 Hz
Temp. Control Range	Amb. ~180° C
Temp. Setting/Display	Digital Control (Controlled by the main controller)

*Double jacketed bath body (Safe to touch)

MODEL	HS-2005S (CE)
Dimension (W × D × H)	720 × 340 × 510 mm
Weight	20 Kg
Power Source	AC 220 V, 2.5 A, 50/60 Hz