

Rotary Evaporator

HS-2005V (CE)

Distinguished Features

6 different types of flasks are attachable (From 50 mL to 3,000 mL)
 Digital speed control
 Digital R.P.M Display/Speed control range : 20~280 R.P.M
 Vertical type condenser
 Digital interior temperature measure and display
 (It catches the temperature changes caused by the interior pressure drops)
 Motorized UP-DOWN control switch
 Easy to utilize standard distillation solution
 Auto timer function
 Autonomous condenser lifting function after a complete shut down of the machine


Included Water Bath HS-3001

(PTFE Coating)

If the amount of water goes down under the certain threshold in the heating bath, the system automatically turns off its power for the safety purpose

Included vacuum controller

*Vacuum range : 0~760 mmHg

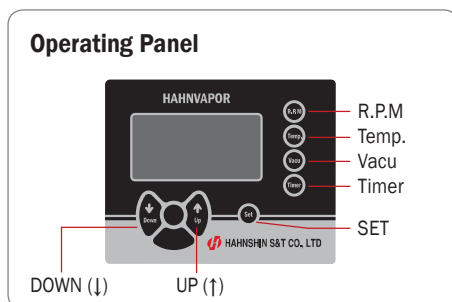

Registration No.

20-0461039/30-0668181

HS-2005V (CE)

Heater	1,000 W
Bath Material	SUS 304
Weight	2.8 Kg
Bath Dimension (Ø × H)	245 × 120 mm (Capacity : 4.3 Liter)
Required Power Source	AC 220 V, 3.5 A, 50/60 Hz
Temp. Control range	Amb. ~180° C
Temp. Setting/Display	Digital Control (Controlled by the main controller)

*Double jacketed bath body (Safe to touch)

MODEL	HS-2005V (CE)
Dimension (W × D × H)	610 × 340 × 730 mm
Weight	21 Kg
Power Source	AC 220 V, 2.5 A, 50/60 Hz